


Journey™ Series – AntiThrust Wedge

Wheelchair Cushion

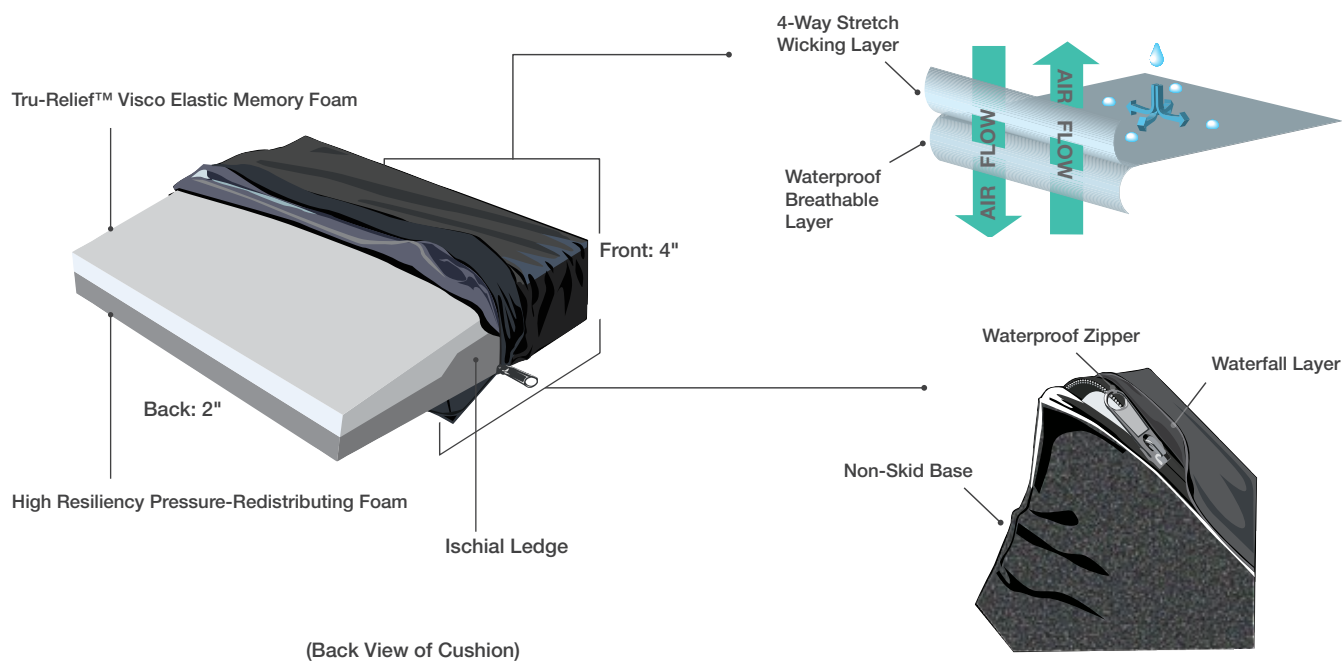

Are the hips to the back of the wheelchair?

If not, consider a:

1. Contoured cushion
2. Shorter length chair (two fingers from edge to popliteal space)
3. Increase hip flexion angle with a wedge or a anti-thrust cushion
4. 45° angle waist belt
5. 90° angle lap belt

JourneyTM Series – AntiThrust Wedge

Wheelchair Cushion

Product Specifications

Wound Risk Level	Medium
Wound Stage	Up to Stage II
Cushion Thickness / Weight Limit	4" front – 2" back / 300lbs

Product Description

The Journey AntiThrust Wedge cushion is designed to offer comfort and positioning. The top layer of Tru-Relief™ visco elastic memory foam conforms to bony prominences. A combination of wedge and AntiThrust design helps to prevent more aggressive forward thrusting and sacral sliding. This combination helps to reduce extensor tone by increasing hip flexion.

Models and Sizes Available (WIDTH x DEPTH x HEIGHT)

Model #	Size
STW4 – 16 x 16	16" x 16" x 4" front to 2" back
STW4 – 16 x 18	16" x 18" x 4" front to 2" back
STW4 – 18 x 16	18" x 16" x 4" front to 2" back
STW4 – 18 x 18	18" x 18" x 4" front to 2" back
STW4 – 20 x 18	20" x 18" x 4" front to 2" back

Comfort • Safety • Mobility

Toll-Free 1(866)330-5336

Keen[®]

Copyright © 2007, Keen Mobility Company. All rights reserved. Rev 070312. All trademarks are property of their respective owner.

Seating 69