


Journey™ Series – Bariatric

Wheelchair Cushion

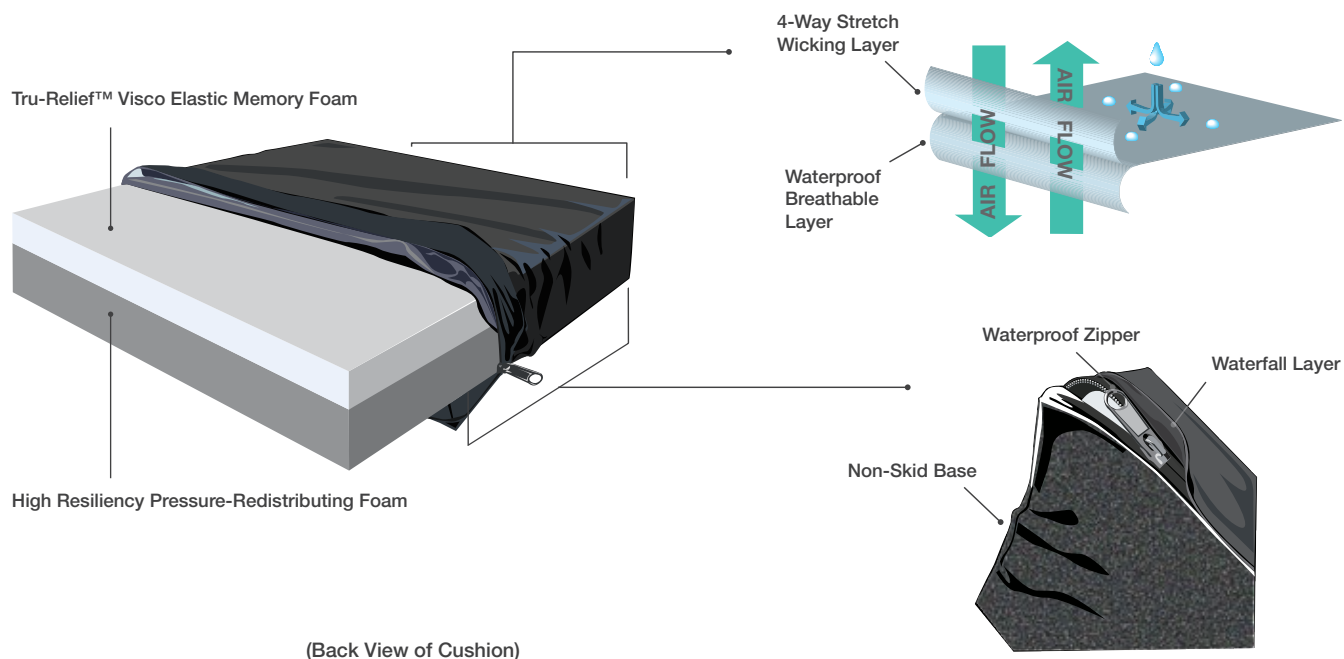

If the seat depth is too short, it may:

1. Increase the risk for an individual to be pushed out of the wheelchair during wheelchair propulsion.
2. Increase pressure over ischial areas because there is less surface area to distribute the weight.

If the seat depth is too long, it may:

1. Restrict circulation.
2. Increase pressure to the back of the knee and calf.
3. Restrict the ability to sit back fully in the seat and receive proper support from the back rest (this results in poor positioning and a “slumped” posture).
4. Inhibit mobility of foot propellers due to lack of clearance at the popliteal space.

