


Journey™ Series – Contour

Wheelchair Cushion

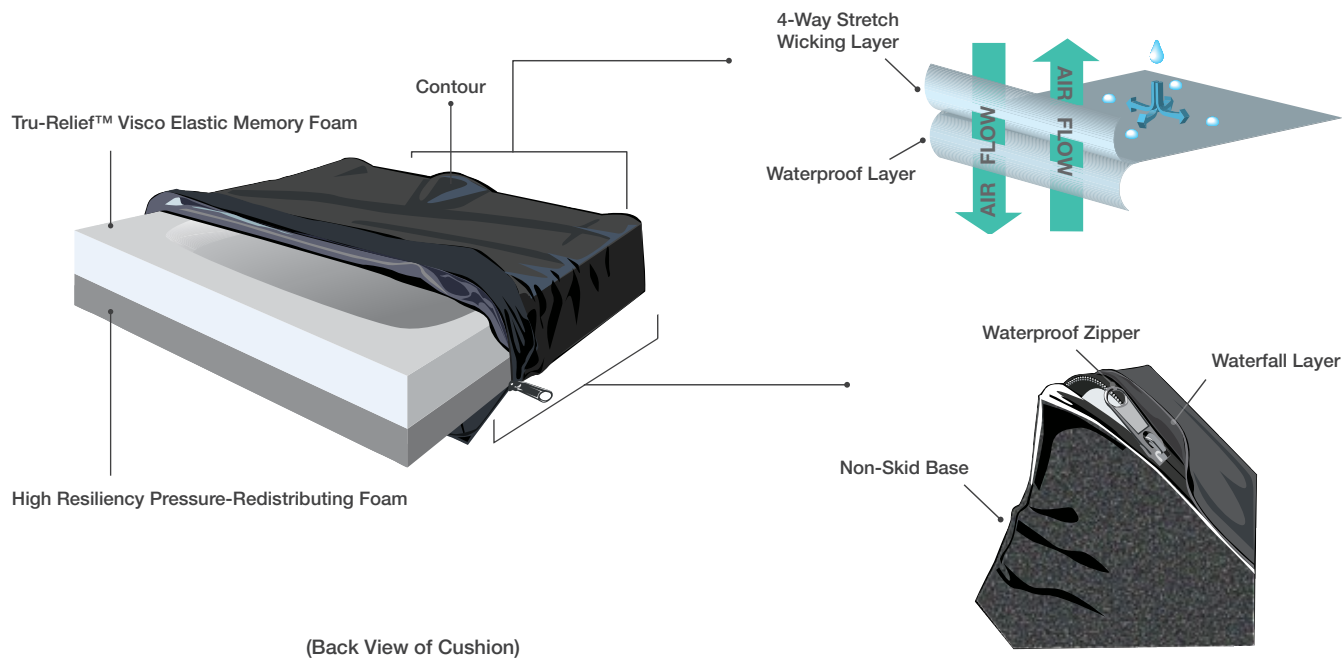

1. Seat width is the first and most important measurement.
2. To ensure proper seat width, measure the hips at the widest part and add 1" on each side.
3. Place both palms between the person's hips and side panels (the hands should contact both hips and side panels without undue pressure or limitation).

Product Specifications
HCPC E2603 (Skin Protection)

Wound Risk Level:	Medium / High	
Wound Stage:	Up to Stage II	
Cushion Thickness / Weight Limit:	2.5" / 250lbs	3.5" / 275lbs

Product Description

The Journey Contour, with more Tru-Relief™ visco elastic memory foam than the Journey Comfort, is designed for individuals at medium to high risk for pressure wounds up to stage II. The slight wedge helps provide mild positioning. The convex bottom is available on the standard 18x16 size to fill in sling seats. The Journey Contour is maintenance-free.

Models and Sizes Available (WEIGHT x DEPTH x HEIGHT)

Contour 2.5"		Contour 3.5"	
Model #	Size	Model #	Size
C2 – 14 x 14	14" x 14" x 2.5"	C3 – 16 x 16	16" x 16" x 3.5"
C2 – 14 x 16	14" x 16" x 2.5"	C3 – 16 x 18	16" x 18" x 3.5"
C2 – 16 x 16	16" x 16" x 2.5"	C3 – 16 x 20	16" x 20" x 3.5"
C2 – 16 x 18	16" x 18" x 2.5"	C3 – 18 x 16	18" x 16" x 3.5"
C2 – 16 x 20	16" x 20" x 2.5"	C3 – 18 x 18	18" x 18" x 3.5"
C2 – 18 x 16	18" x 16" x 2.5"	C3 – 18 x 20	18" x 20" x 3.5"
C2 – 18 x 18	18" x 18" x 2.5"	C3 – 20 x 18	20" x 18" x 3.5"
C2 – 18 x 20	18" x 20" x 2.5"	C3 – 20 x 20	20" x 20" x 3.5"
C2 – 20 x 18	20" x 18" x 2.5"		
C2 – 20 x 20	20" x 20" x 2.5"		
convex bottom	convex bottom		
CC2 – 18 x 16	18" x 16" x 2.5"		

Comfort • Safety • Mobility
Toll-Free 1(866)330-5336


Copyright © 2007, Keen Mobility Company. All rights reserved. Rev 070312. All trademarks are property of their respective owner.

Seating 51