


# Journey™ Series – Contour Plus

Wheelchair Cushion

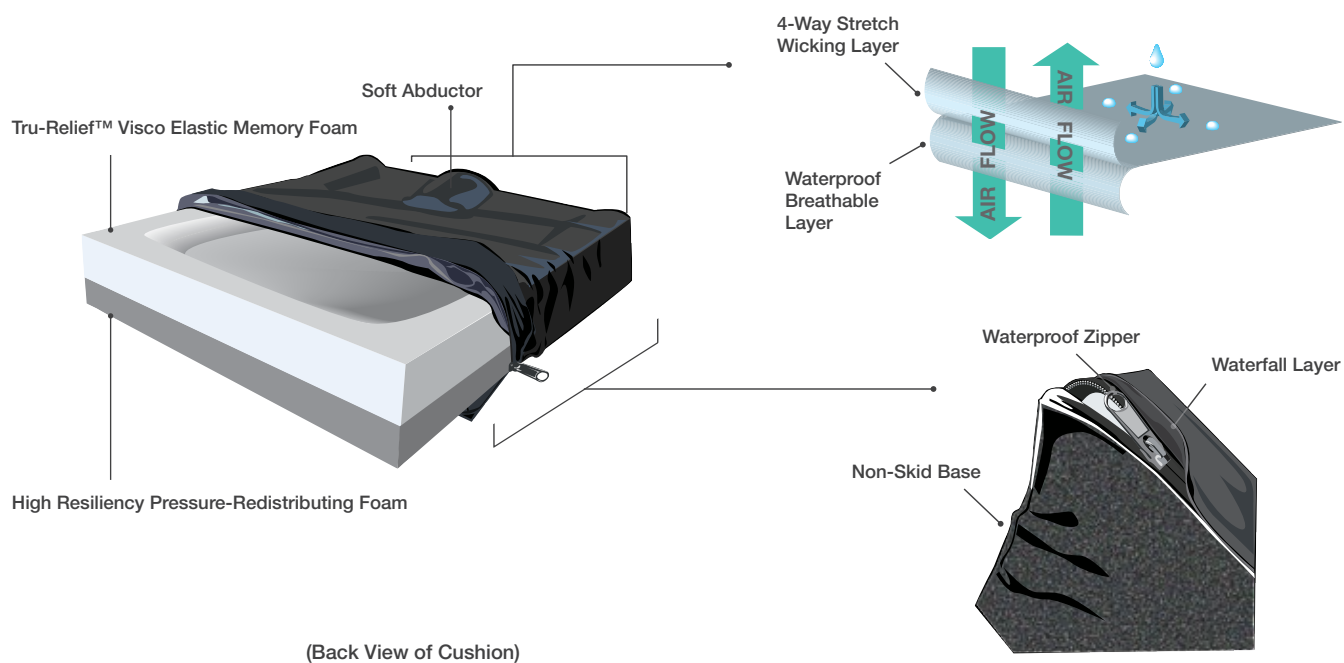

If seat width is too narrow:

- An individual has the potential for developing pressure areas due to increased pressure on the trochanters and ischium.
- There may be transfer difficulties since mobility is limited.

# Journey<sup>TM</sup>Series – Contour Plus

Wheelchair Cushion

## Product Specifications HCPC E2607 (Skin Protection / Positioning) & E2608 (22" or greater)

Wound Risk Level:	Medium / High	
Wound Stage:	Up to Stage III	
Cushion Thickness / Weight Limit:	2.5" / 250lbs	3.5" / 275lbs

## Product Description

The Journey Contour Plus cushion combines Tru-Relief™ visco elastic memory foam with moderate positioning. Designed to encourage proper pelvic alignment and a neutral spine, the deep contour, abductor and slight wedge helps prevent the pelvis from posterior pelvic tilt and sacral sliding. The maintenance-free Contour Plus Cushion is a great solution for individuals with Stage II & III wounds or a history of wounds.

## Models and Sizes Available (WEIGHT x DEPTH x HEIGHT)

Contour Plus 2.5"		Contour Plus 3.5"	
Model #	Size	Model #	Size
CP2 – 16 x 16	16" x 16" x 2.5"	CP3 – 16 x 16	16" x 16" x 3.5"
CP2 – 16 x 18	16" x 18" x 2.5"	CP3 – 16 x 18	16" x 18" x 3.5"
CP2 – 18 x 16	18" x 16" x 2.5"	CP3 – 18 x 16	18" x 16" x 3.5"
CP2 – 18 x 18	18" x 18" x 2.5"	CP3 – 18 x 18	18" x 18" x 3.5"
CP2 – 18 x 20	18" x 20" x 2.5"	CP3 – 18 x 20	18" x 20" x 3.5"
CP2 – 20 x 16	20" x 16" x 2.5"	CP3 – 20 x 16	20" x 16" x 3.5"
CP2 – 20 x 18	20" x 18" x 2.5"	CP3 – 20 x 18	20" x 18" x 3.5"
CP2 – 20 x 20	20" x 20" x 2.5"	CP3 – 20 x 20	20" x 20" x 3.5"
CP2 – 22 x 16	22" x 16" x 2.5"	CP3 – 22 x 16	22" x 16" x 3.5"
CP2 – 22 x 18	22" x 18" x 2.5"	CP3 – 22 x 18	22" x 18" x 3.5"

**Comfort • Safety • Mobility**

**Toll-Free 1(866)330-5336**

**Keen**<sup>®</sup>

Copyright © 2007, Keen Mobility Company. All rights reserved. Rev 070312. All trademarks are property of their respective owner.

**Seating 53**